INFINITY5-5

5.1 Megapixel Global Shutter CMOS Microscope Camera

High performance microscope camera for a wide range of applications – with dual output to HDMI and USB 3



INFINITY5-5

Lumenera's INFINITY5-5 is a high quality 5 megapixel microscopy camera that provides high speeds at high resolution. The INFINITY5-5 is based on the Sony® Pregius™ global shutter CMOS sensor that rivals sCMOS technology. With fast focusing at high frame rates of up to 75 frames per second, the INFINITY5-5 performs in a wide range of applications.

The INFINITY5-5's 2/3 inch sensor format accommodates a 2464 x 2056 resolution with 3.45 micron pixels. The HDMI output allows the INFINITY5-5 to connect to both a computer and HDMI TV or monitor for applications where knowledge sharing is critical.

Simultaneous Live Viewing on Computer and HDMI Monitor

Operate the camera by directly connecting to an HDMI monitor, or connect the camera to a PC or Mac via USB 3 for image preview and capture. Use the soft-touch response buttons to control the camera with or without connection to a computer.

Ready to Perform


The INFINITY5-5 comes ready to perform with cables (USB and HDMI), power supply, and INFINITY CAPTURE software for routine imaging needs, including camera setting adjustments, live preview, image capture, and video clips. The INFINITY5-5 also provides third-party interoperability with MetaMorph and Micro-Manager microscopy software, MatLab systems engineering software, LabView analysis software, and TWAIN and DirectShow interfaces.


Accurate Color

Lumenera's specialized algorithms ensure accurate color reproduction that results in a sample preview on a monitor that matches the view in the microscope eyepiece.

Guaranteed Quality

Lumenera provides an industry-leading, full four-year year warranty on all microscopy cameras.

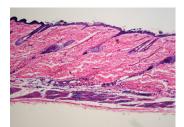

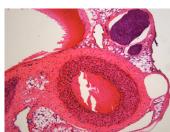




Features

INFINT.


Lumenera

- Sony Pregius global shutter sensor technology
- Color or monochrome IMX250 Global Shutter CMOS sensor with 2/3" optical format providing a resolution of 2464 x 2056 using 3.45 x 3.45µm pixels
- Fast frame rates at ~75 fps at full resolution in 8-bit
- Simultaneous output to USB 3 and HDMI compatible monitor
- Intuitive on-camera buttons for power, white balance and image capture
- Selectable 8 or 12-bit pixel data
- High-speed USB 3.1 Gen 1 interface for fast image delivery and simplified connectivity. USB 2.0 supported.
- Software compatible with Windows 10, Windows 8.1, Windows 7, MAC OS X 10.13, 32 and 64-bit operating systems
- Support for third-party capture and analysis applications: MetaMorph and Micro-Manager
- Recommended coupler: 0.6x 0.7x

Color Quantum Efficiency

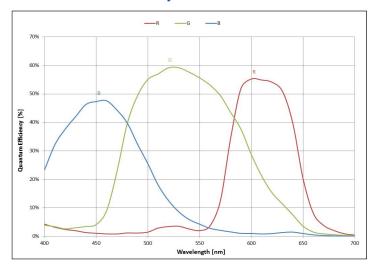

Monochrome Quantum Efficiency

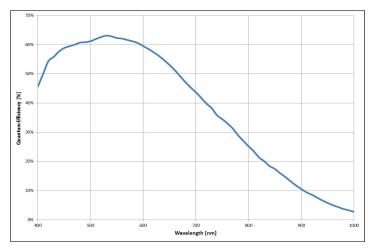

Ordering Information	
INFINITY5-5C	5.1 MP CMOS Color USB 3.1 Gen 1 Camera
INFINITY5-5M	5.1 MP CMOS Mono USB 3.1 Gen 1 Camera
LuIAP-2	INFINITY Advanced Features Pack 2: Includes USB Key for extra INFINITY ANALYZE license + Advanced Features Module, 5 year total warranty, 1 advance product replacement
La050315	Transformer, 5VDC, 15W, 3A, International
La2000PAFL	GPIO cable with leads
La2030HD	3m HDMI Cable

Sensor Specifications		
Image Sensor	SONY IMX250 CMOS color or monochrome	
Optical Format	2/3"	
Imager Size	Diagonal 11.1 mm	
Pixel Size	3.45 x 3.45 μm	
Resolution	2464 x 2056 pixels	
Region of Interest Control	Yes	
Camera Specifications		
Max Frame Rate	1080P60 (~75 fps at full resolution) in 8-bit	
Bit Depth	8 or 12-bit	
Binning Modes	2 x 2 for mono	
Exposure Control	Manual and automatic control	
Exposure Range	42µs to 59.5m (snapshot); 14µs to 14s (video)	
Gain Control	Manual and automatic control	
Gain Range	~1 to 16 x	
White Balance	Manual and automatic control	
Trigger Modes	Hardware and software triggerable	
Camera Characteristics	Hardware and software inggerable	
	Mono: 4.9 DN/(nJ/cm2), Color: 4.3 DN/(nJ/cm2) (Global	
Sensitivity (8-bit)	and channel gains at unity)	
Dynamic Range	~72 dB	
Full Well Depth	~10,500 e-	
Quantum Efficiency	59%@ 530nm (color), 63%@ 530nm (mono)	
Read Noise	~2.30 e-	
Dark Current	1.0 e-/s (at 22 °C ambient, 35 °C internal camera)	
Mechanical Specifications		
Data Interface	USB 3 connector, HDMI output	
General Purpose I/O	Locking Hirose MXR-8R-8SA(71)	
Lens Mount	Adjustable C-Mount	
Dimensions	3.85 x 2.75 X 2.00 inch	
Mass	330 g 0 to 50 °C	
Operating Temperature	-30 to 70 °C	
Storage Temperature Operating Humidity	5 to 95 %, non-condensing	
Shock / Vibration	50 G Shock, 5 G (2-200Hz) Vibration	
Camera Software		
	Windows 10, Windows 8.1, Windows 7, MAC OS X	
Operating Systems	10.13, 32 and 64-bit operating systems	
Third Party Software Interfaces	MicroManager, MetaMorph, MatLab, LabView, Direct- Show, TWAIN	
Power and Emissions		
Power Consumption	~4W	
Power Requirement	External 5 V DC, 1.2 A, power supply (included)	
Emissions Compliances	FCC and CE Certified	
Hazardous Materials	RoHS, WEEE Compliant	
Warranty	Four (4) years	
Included In The Box		
INFINITY <i>5-</i> 5	5.1 MP digital camera with 3m USB 3 cable	
LuINFLTSW-CD	CD with INFINITY ANALYZE and CAPTURE software, TWAIN driver and documentation	
La050315	Transformer, 5VDC, 15W, 3A, International	
La2030HD	3m HDMI Cable	



Mechanical Drawing

